




SOLVE EVERY CASE.
WITH 100% CERAMIC.

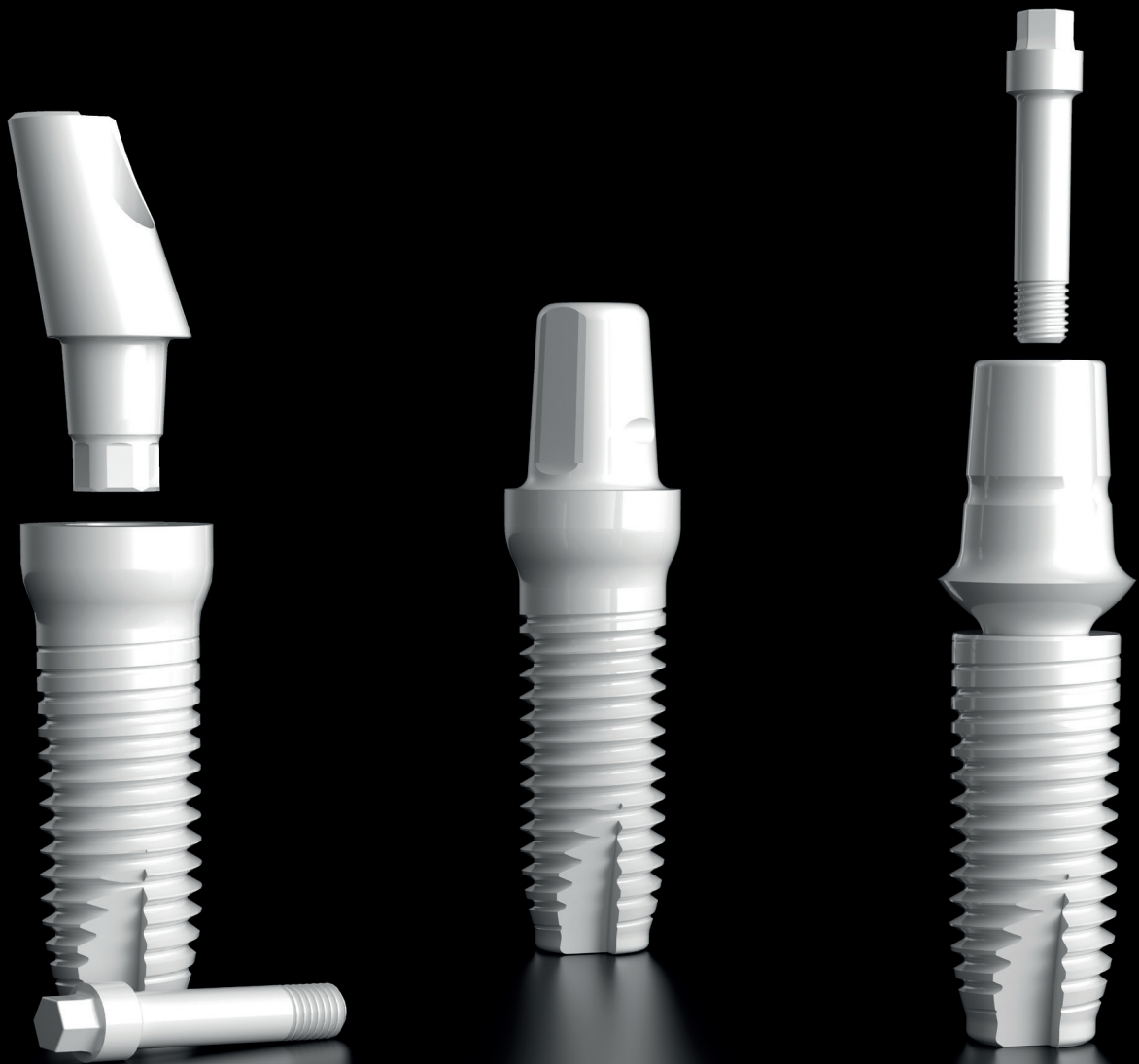

Global leader in ceramic implants. Since 2004.

Z-SYSTEMS is the global leader in ceramic implants.

Tens of thousands of successfully placed implants and over 15 years' experience have brought us significant advances in manufacturing, quality and know-how. We value clinical documentation, and are certified as per ISO 13485.

Z-SYSTEMS products meet the world's highest certification standards.

We obtained the first CE certification for one-piece ceramic implants in 2004. All our ceramic implants are manufactured with our safe Zirkolith® process. In 2007, we were the first ceramic implant manufacturer to meet the requirements to enter the demanding US market. Four years later, we continued this success with our two-piece ceramic implant.

Since 2016, we are the first and only manufacturer of a screw-retained bone-level implant made of ceramic with platform switching.

Then, in 2018 with the Z5-BL and in 2020 with the Z5-TL, Z-SYSTEMS has virtually revolutionized ceramic implantology. The world's first two-piece screw-retained implants entirely without plastic or metal components enable users to perform workflows and indications that were previously only conceivable with titanium implants.


Trusted pioneering advances.

Constant innovation, based on solid scientific experimentation, our unique product range and personal advice – your satisfaction is our goal. The road to reliable, safe material and design is a long one. We follow significant quality demands and continue successive developments. Safety without compromise is our highest motto, because we are convinced that only sustainable and safe implants ensure your success.

We are pleased to present you our complete product range on the following pages, and to win you over with our safe ceramic implants.

Your Z-SYSTEMS Team

- ✓ Survival rate comparable to that of titanium implants
- ✓ Secure osseointegration thanks to the hydrophilic SLM® surface
- ✓ Implant and abutment can be prepped
- ✓ Excellent biocompatibility
- ✓ Less plaque, better gum attachment
- ✓ Not electrically conductive, no galvanic elements, no currents
- ✓ Sustained aesthetic results

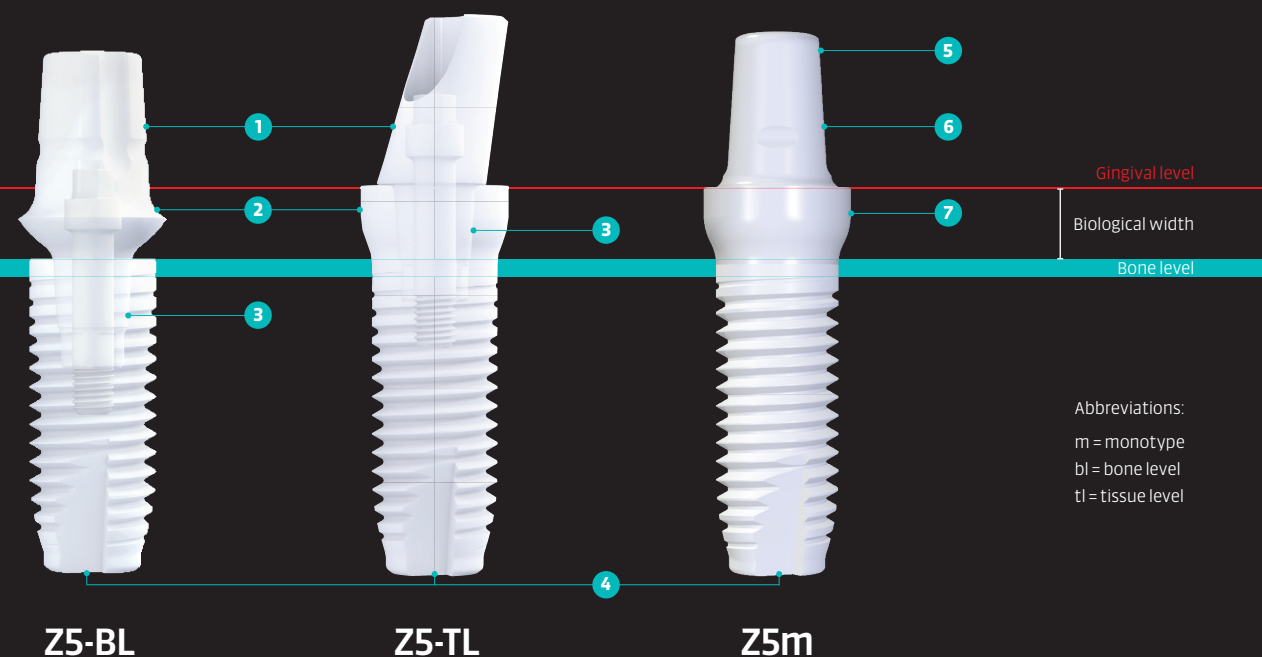

The most advanced & versatile ceramic dental implants on the market today

- 1 100% ceramic, screw-retained, two-piece Bone/Tissue Level Implant: Surgical procedure analogous to two-piece titanium implantology.
- 2 Straight and angled abutments in different abutment and gingival heights are available for aesthetically superior restorations.
- 3 The extraordinarily stable implant-abutment connection plus outstanding break resistance ensure the highest level of safety.
- 4 The conical, self-cutting threads ensure primary stability.
- 5 Z5m can be prepped.* Zirkolith® material's high stability makes it possible to prep the implant to defined preparation limits. This unique ability in implantology means a great deal of freedom for the user for individual adjustment to existing anatomical circumstances.
- 6 Our abutment connection geometries are integrated and accessible within the most widely used CAD libraries, including Exocad, 3Shape, Swissmeda/SMOP and SiCat. This ensures seamless digital workflow integration for dental professionals, ensuring efficient planning and precise design.
- 7 This unique material makes it possible to prep the shoulder on a custom basis to defined preparation limits.* The shoulder design makes it possible to produce the optimal soft tissue shape for a safe and aesthetic solution.
- 8 The optimal combination of the conical basic form with dynamic, active and self-cutting threads achieves high primary stability even in soft D3 bone.

All Z-SYSTEMS implants are provided with a hydrophilic surface! The laser-structured SLM®-surface has been shown to promote osseointegration due to its demonstrable increased roughness.

* Z-SYSTEMS has had the safety of manually-prepped Z-SYSTEMS implants tested by the renowned Fraunhofer Institute for Material Engineering. The results are clear. No significant impairment was found in the prepped implants. Manually prepped implants do not have a higher risk of breakage.

Z5-BL The Bone Level implant

The world's first 100% ceramic, two-piece, screw-retained implant-abutment connection

Z5-BL 40-10


> Implant 4 mm × 10 mm

100% ceramic

- 100% ceramic, screw-retained, two-piece bone level implant
- Surgical procedure analogous to two-piece titanium implantology
- Closed healing, augmentation and simultaneous sinus lift is possible in suitable cases

Experience from over 80,000 implants placed

- Established enossal design, based on tissue-level predecessors, ensures transferable results
- Self-tapping thread with bone chip reservoir reduces insertion torque
- Excellent osseointegration through Z-SYSTEMS's proven, patented proprietary SLM® laser surface technology and gentle plasma sterilization

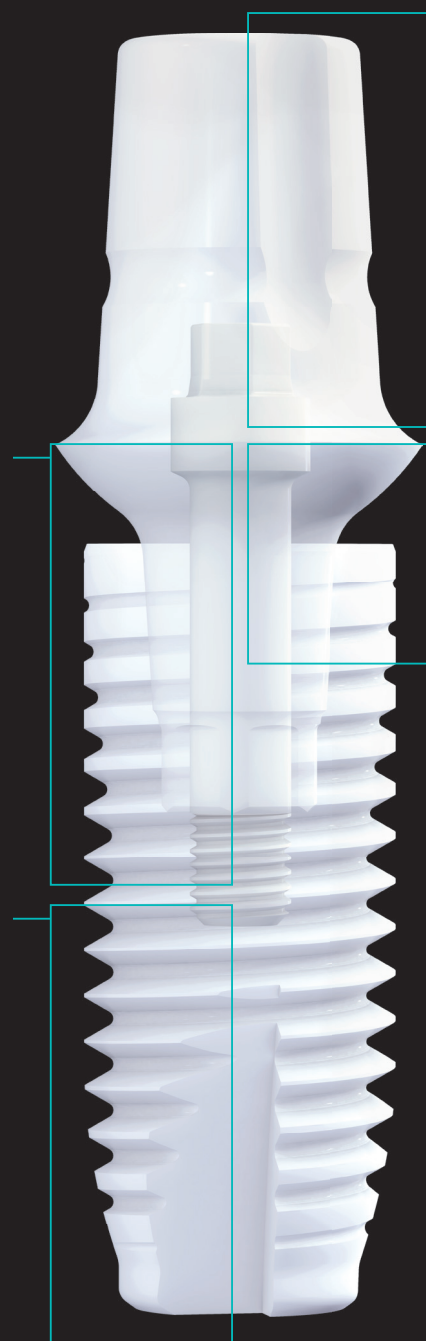


Aesthetics for professional care

- Wide range of straight and angled abutments for optimal esthetic results and flexible functionality
- A wide range of options for designing the exit profile

A unique, strong connection

- The conical screw connection with internal thread for micro-gap-free connection and prevention of a pump effect.
- NEW:** ceramic abutment screw for 100% all-ceramic restoration
- A titanium screw is also available as an option
- Platform switch for a more stable and enlarged mucosal cuff and for reduced loading of the crestal bone

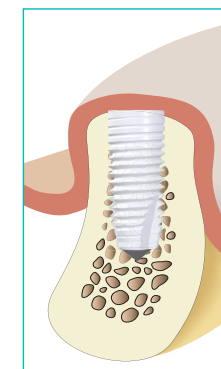

Z5-BL

The new screw-retained Bone Level Implant

	D 3.6 mm	D 4.0 mm	D 5.0 mm
8	Z5-BL-3608	Z5-BL-4008	Z5-BL-5008
10	Z5-BL-3610	Z5-BL-4010	Z5-BL-5010
12	Z5-BL-3612	Z5-BL-4012	Z5-BL-5012

■ = length in mm

Abutments (incl. screw of your choice and single-use screwdriver)




BL-A0015, BL-A0025, BL-A1515, BL-A1525, BL-AN1515, BL-BB1545, BL-BB2545, BL-CB1545, BL-CB2545, BL-CB1555, BL-CB2555, BL-CB2565




BL-LC0030, BL-LC0040, BL-LC1530, BL-LC1540, BL-ZB1538, BL-ZB1545, BL-ZB2545, BL-TA0030, BL-TA1538

Healing cap

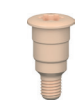

BL-HC

Gingiva formers

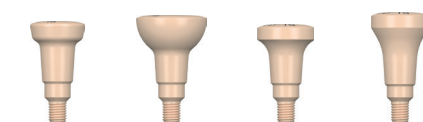


BL-GF1538, BL-GF3055, BL-GF1545, BL-GF2545



BL-GF1555, BL-GF2555, BL-GF2565

Impression tools

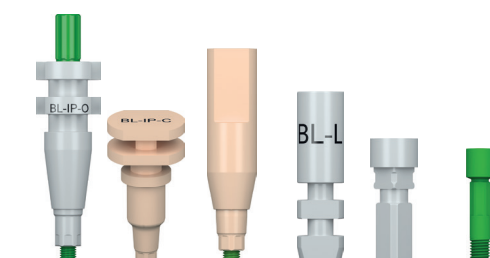

BL-IP-O, BL-IP-C, BL-SB-36*, BL-L, BL-L-3D*, BL-OSL

*Digital workflow currently only available with exocad

Screws



BL-OSC-H, BL-OST-H KI589, ZT-HA-9, ZT-RA10-9, ZT-RA20-9, TP-SNS, BL-SD-ST, BL-SD-LT, SD-B5-S, SD-B5-L, SD-T6-S, SD-T6-L, BL-CD

Bone taps



T360-3, T400-3, T500-3

Z5-TL

The Tissue Level implant

The world's first 100% ceramic, two-piece, screw-retained implant-abutment connection

Z5-TL 40-10

> Implant 4 mm × 10 mm

100% ceramic

- 100% ceramic, screw-retained, two-piece Tissue Level Implant
- Surgical procedure analogous to two-piece titanium implantology
- Optimized soft tissue management thanks to patented shoulder design

Experience from over 80,000 implants placed



- Established enossal design, based on tissue-level predecessors, ensures transferable results
- Self-tapping thread with bone chip reservoir reduces insertion torque
- Excellent osseointegration through Z-SYSTEMS's proven, patented proprietary SLM® laser surface technology and gentle plasma sterilization

Aesthetics for professional care

- Wide range of abutments for optimal esthetic results and flexible functionality
- A wide range of options for designing the exit profile

A unique, strong connection

- The conical screw connection with internal thread for micro-gap-free connection and prevention of a pump effect.
- NEW:** ceramic abutment screw for 100% all-ceramic restoration
- A titanium screw is also available as an option

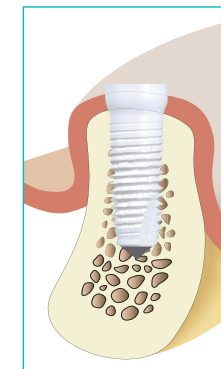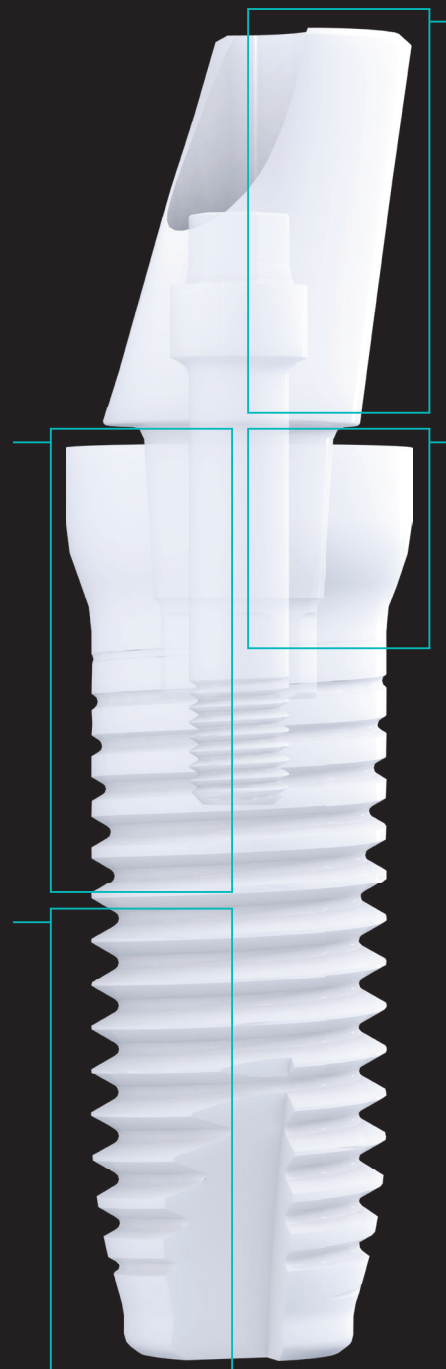

Z5-TL

The new screw retained Tissue Level Implant


	D 4.0 mm	D 5.0 mm
8	Z5-TL-4008	Z5-TL-5008
10	Z5-TL-4010	Z5-TL-5010
12	Z5-TL-4012	Z5-TL-5012

■ = length in mm

Abutments (Incl. screw of your choice and single-use screwdriver)




Healing cap

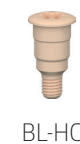


BL-HC

Gingiva formers

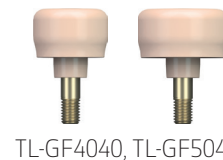

TL-GF4040, TL-GF5040

Impression tools

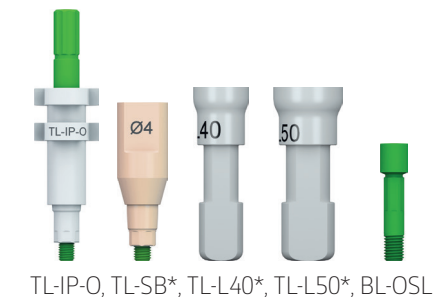

TL-IP-O, TL-SB*, TL-L40*, TL-L50*, BL-OSL

*Digital workflow currently only available with exocad

Screws



BL-OSC-H, BL-OST-H

Instruments



KI589, ZT-HA-9, ZT-RA10-9, ZT-RA20-9, TP-SNS, BL-SD-ST, BL-SD-LT, SD-B5-S, SD-B5-L, SD-T6-S, SD-T6-L, BL-CD

Bone taps



T400-3, T500-3

Z5m Monotype implant


The trusted implant

Z5m-40-10

> Implant 4 mm × 10 mm


Experience from over 60,000 implants placed

- Excellent osseointegration through Z-Systems's proven, patented proprietary SLM® laser surface technology and gentle plasma sterilization
- Self-tapping thread with bone chip reservoir reduces insertion torque

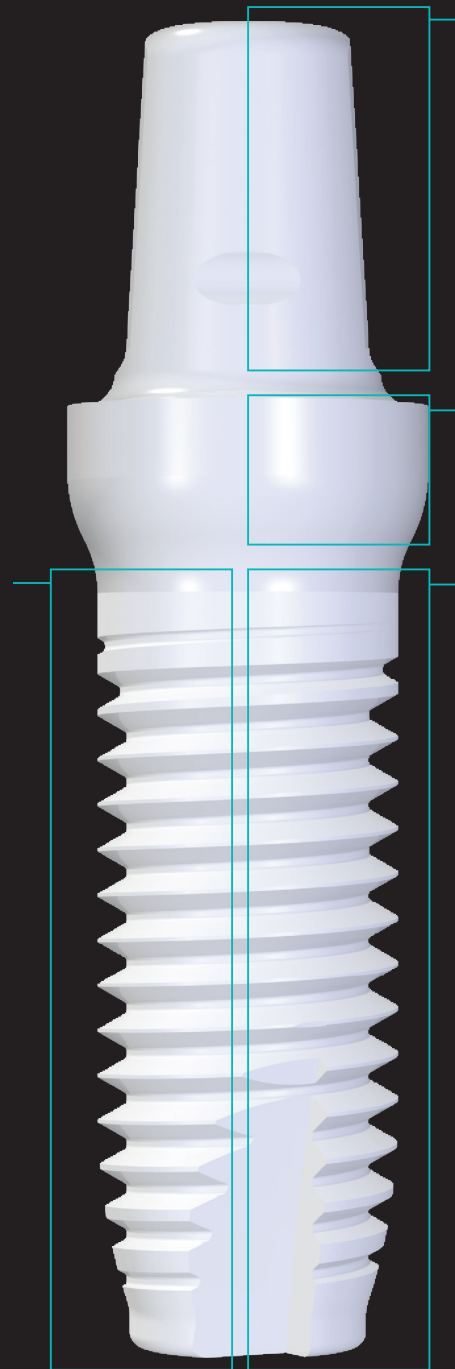

Can be prepped directly in situ: Abutment and implant shoulder

- individual aesthetics
- Fast and cost-effective restorations – like a natural tooth
- Perfect individual aesthetics due to preppable preparation margin.
- Conically shaped thread with broadened core diameter in the upper thread area ensures optimal mechanical strength and better primary stability

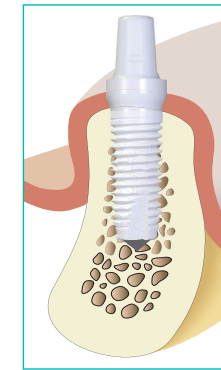

Z5m

The trusted implant
for a variety of indications

	Ø 3,6 mm	Ø 4,0 mm	Ø 5,0 mm
8		Z5m-40-08	Z5m-50-08
10	Z5m-36-10	Z5m-40-10	Z5m-50-10
12	Z5m-36-12	Z5m-40-12	Z5m-50-12

■ = Length in mm


Gingiva formers



Impression caps



Laboratory analogues

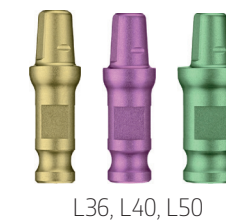


Bone taps




Instruments



Drilling protocol Z5-BL

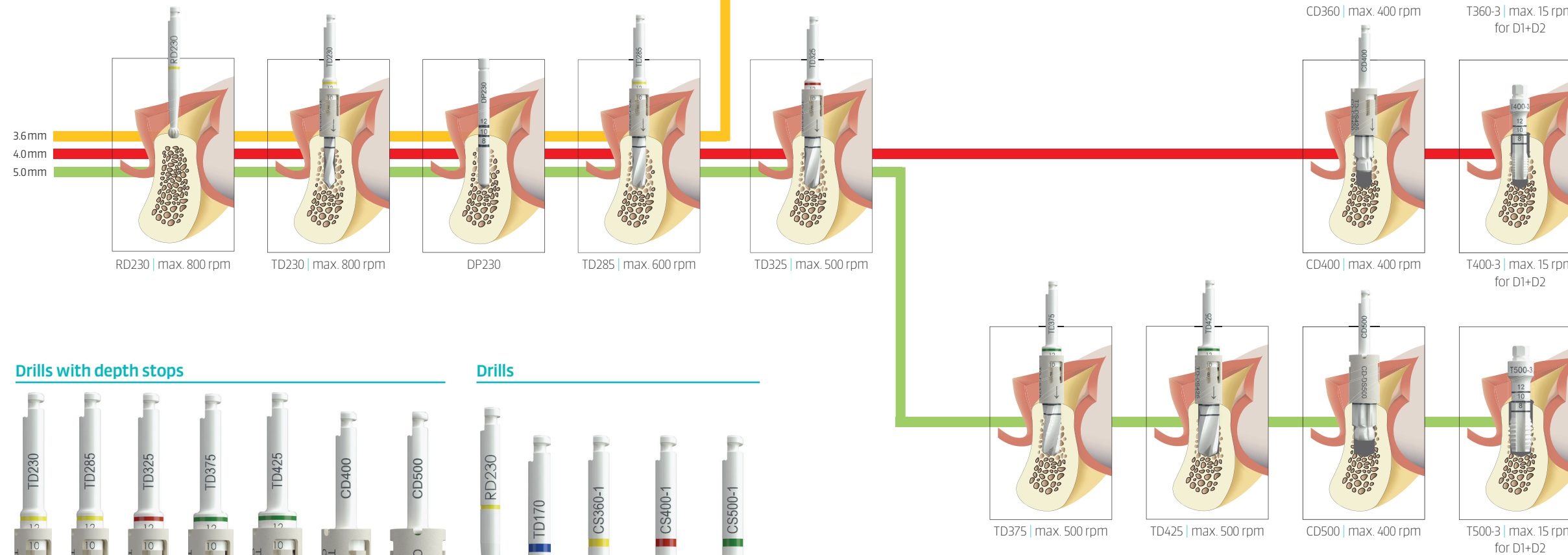

Drills with depth stops

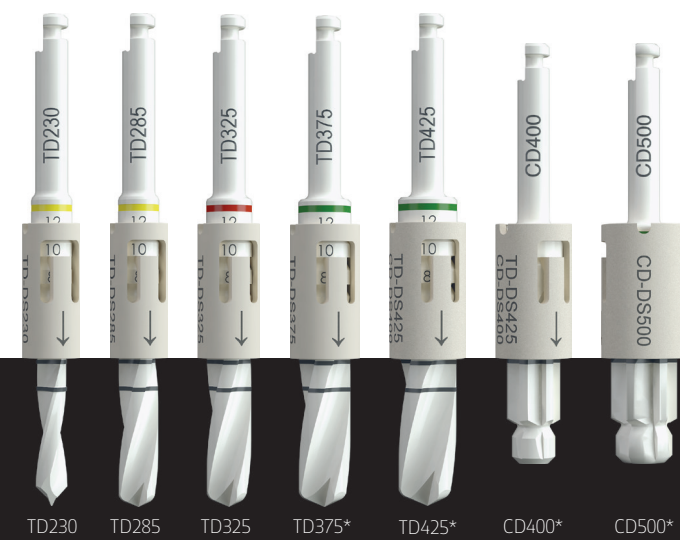


Drills




*Used only in certain site preparations. See user manual of your selected implant type.


Drilling protocol Z5-TL

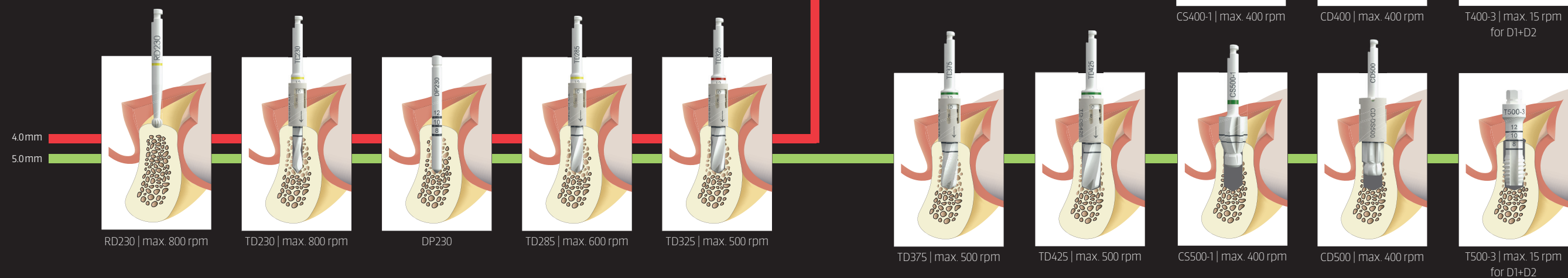


Drilling protocol Z5m

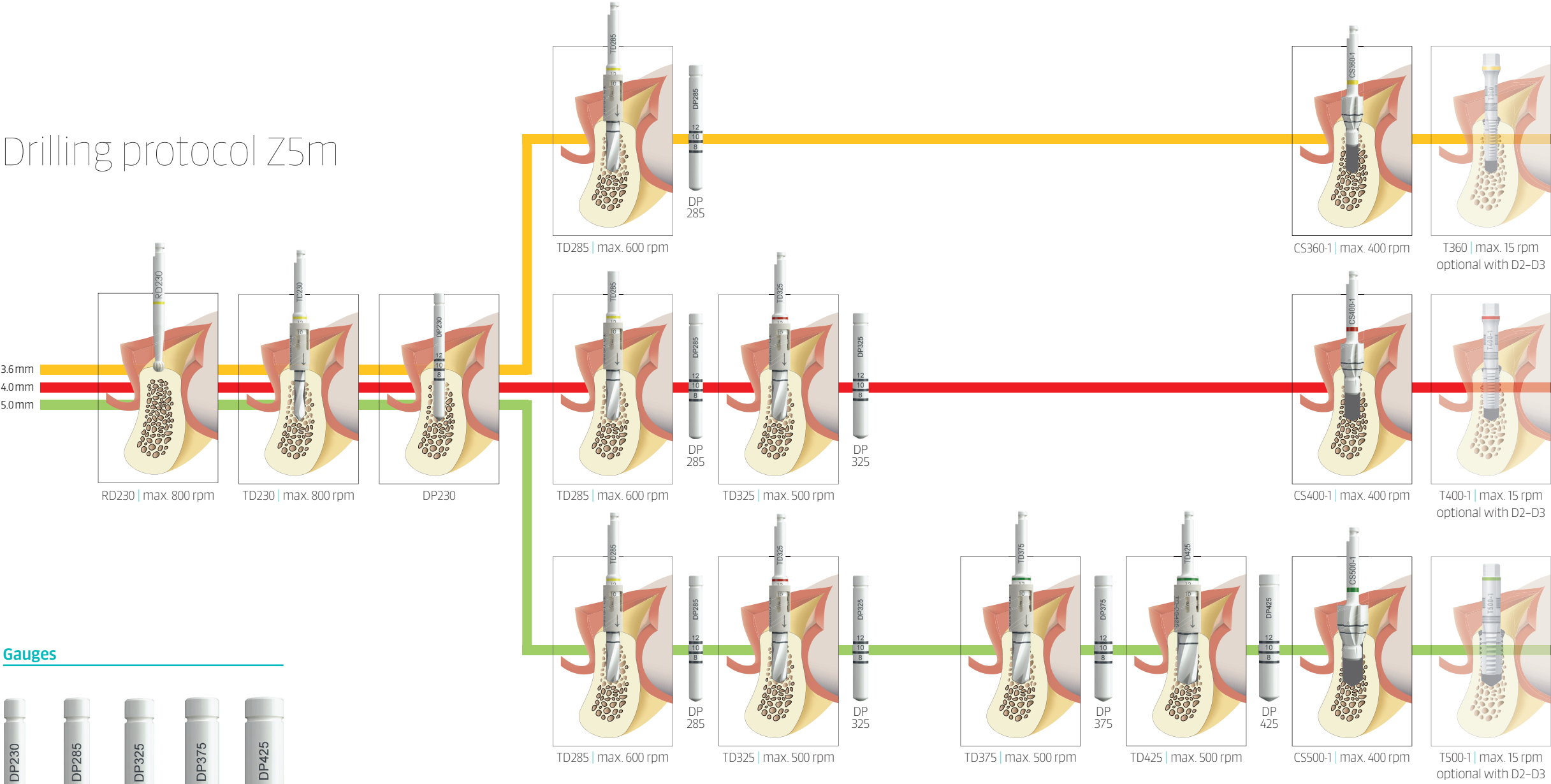


Gauges

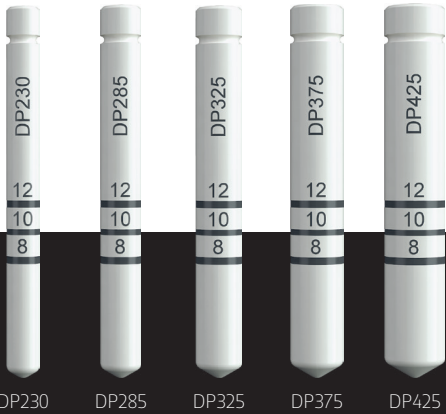

Drilling depth for Z5m implants according to the implant length: 8/10/12 mm

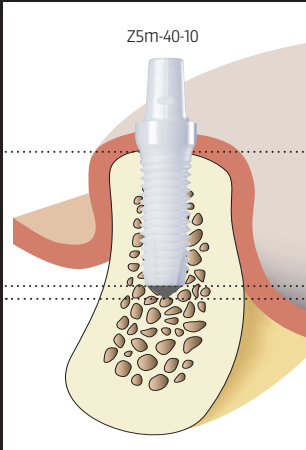

Insertion Depth

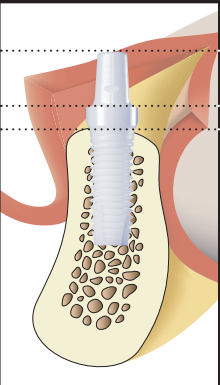

Abutment 5mm
Biolog. width 2.5mm
bone level

What our Doctors are saying

“

Having used Z-SYSTEMS all-zirconia implants for more than nine years, I have seen this company's commitment to continuously improving their products. They are currently introducing their sixth generation of zirconia implants to the U.S., when no other company is past their first generation zirconia implant in the U.S. market. In fact, most other companies haven't even been in the market for a year.

I look fondly at our 8-year results. The implants, the restorations, and the soft tissues look as beautiful as the day the implants were restored. There is no better motivation to go back again and again to Z-SYSTEMS products for each new case!



Dr Ted Fields DDS PhD
Dallas, USA

”

“

Z-SYSTEMS has opened up a whole new world of patient solutions in ceramic implants, with their full line of a 1-piece, 2-piece & Tapered Implants, it's no wonder they lead the market with the most trusted and tested full line of Ceramic Dental Implants. A proven alternative to Titanium Implants: finally making metal-free dentistry a reality



Dr John Roberts, B.Ch.D
Huddersfield, GBR

”

“

Z-SYSTEMS ceramic implants offer a revolutionary, innovative and metal-free alternative to traditional titanium implants and set new standards in dentistry. Their outstanding biocompatibility promotes optimal bonding to the gums, virtually eliminating inflammation. The white color of the implants guarantees natural aesthetics, and ceramic also prevents allergies and intolerances, making them an optimal choice for patients and dentists who value health and aesthetics. Z-SYSTEMS implants are a future-proof, aesthetic and reliable solution in modern implantology.



Dr. med. dent. Judith Th. Aufenanger MSC., MSC.
Küsnacht, CH

”

“

What I like most is that with Z-SYSTEMS' broad product range, actually any clinical situation can be mastered. While one-piece implants still remain a safe value for certain cases, the two-piece screw-retained implant types now available open up even more possibilities. And since many of my patients have been wearing their implants for over 15 years, I know that with Z-SYSTEMS, long-term stability is also a given.



Dr. Ralf Lüttmann
Eckernförde, GER

”

“

Z-SYSTEMS is my first choice in Aesthetic dentistry. For patients who think holistically and care about a harmonious interaction with their bodies, these ceramic implants are simply second to none. Metals tend to trigger reactions in the body, but with the exceptional biocompatibility of Z-SYSTEMS ceramics, you can be sure that your dentures will be optimally tolerated. With Z-SYSTEMS, you not only get aesthetically pleasing results, but also the peace of mind that your health is the main focus.



Kian Dilmaghani M.Sc.
Basel, CH

”

“

Z-SYSTEMS is a true game changer in the implant arena. It is clear why Z-SYSTEMS leads the market with the most trusted ceramic implants in the world. Having used ceramic implants for the past 6 years, I was very excited to start using Z-SYSTEMS products in the last year.

Z-SYSTEMS ongoing commitment to quality & continual innovations provide practitioners with more flexibility and options than any other ceramic implant company, and allows me to treat the specific needs of each and every patient. History matters when it comes to quality treatments!



Dr Al Manesh, Periodontist
Mission Viejo, USA

”





Headquarters

Z-Systems AG | Werkhofstrasse 5 | CH-4702 Oensingen | Tel. +41 62 388 69 69 | support@zsystems.com

International Distribution

Z-Systems GmbH | Wilhelmstraße 6 | D-78120 Furtwangen | Tel. +49 7723 9144-899 | support@zsystems.com

www.zsystems.com